

SmartSite® Needle-Free Access Systems



Simple, effective and proven systems for your entire continuum of care

SmartSite® is one of the most widely used needle-free access systems in the world.

- Smooth top can be effectively swabbed with standard, acceptable disinfection practice
- Straight-through fluid path helps improve the ability to flush
- Zero dead space helps prevent issues caused by blood clotting, difficulty in flushing fluid pathway or interactions with incompatible drugs

SmartSite® Needle-Free Systems

Standardization throughout the hospital

The SmartSite® valve allows you to standardize on one needle-free technology hospital-wide:

- Metal-free
- DEHP-free, lipid-resistant and highly chemical compatible
- Latex-free

Smooth, swabbable surface


The SmartSite® valve has a smooth top that is easy to swab. IV connectors with a smooth design can be effectively swabbed with a standard, acceptable disinfection practice.

Straight-through fluid path


The SmartSite® valve has a straight-through fluid path (approx. 0.1mL) equivalent to an 18g needle/catheter, which helps reduce the risk of hemolysis.

Securement seal

When the SmartSite® valve is activated, there is an airtight seal between the male Luer and the top of the blue piston of the valve. This helps ensure there is no air entrapment or fluid leakage when the valve is activated.

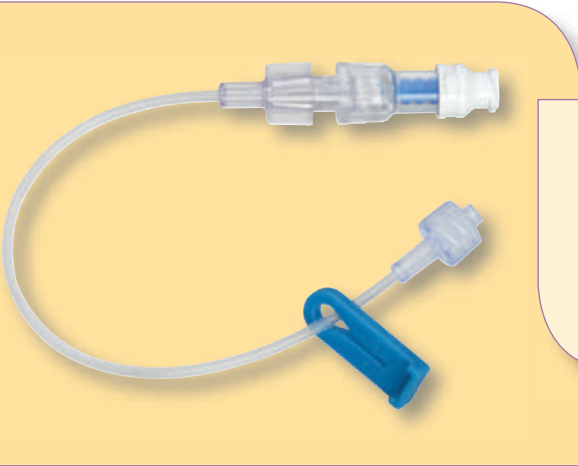

The valve's ergonomic design makes it intuitive and easy to use.



Smiths Medical offers several options of needle-free access systems with SmartSite® valves.

Contact your local Smiths Medical sales representative or customer service for additional information.

Available configurations

- T-connector extension set
- Filter extension set
- Gravity administration set
- Bifurcated extension sets
- Stopcock extension sets
- Manifold extension sets
- Trifurcated extension sets

Smiths Medical ASD, Inc.
St. Paul, MN 55112, USA
Phone: 1-214-618-0218
Toll-Free USA: 1-800-258-5361
www.smiths-medical.com

Find your local
contact information at:
[www.smiths-medical.com/
customer-support](http://www.smiths-medical.com/customer-support)

The Smiths Medical and Medex design marks are trademarks of Smiths Medical. The symbol © indicates the trademark is registered in the U.S. Patent and Trademark Office and certain other countries. All other names and marks mentioned are the trademarks or service marks of their respective owners. ©2012 Smiths Medical. All rights reserved. 9/12 IN193858EN-0912

smiths medical